


## EMERGENCY INLET POINT PANEL FOR CENTRAL STATIONS

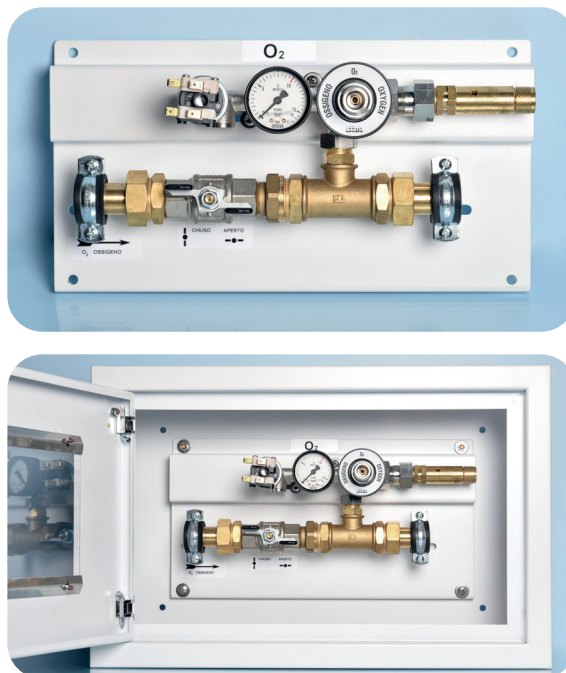

**THE EMERGENCY INLET POINT PANEL** allows to by-pass the station and to keep supplying the network by means a supplying point in case of gas central station panel board maintenance. The unit integrates a safety valve and a pressure-switch for the connection to the alarm unit.

The unit is constituted by a varnished steel frame having built-in:

- .1 sectioning ball valve complete of 3 pieces brass fitting
- .1 primary network pressure-switch or pressure trasducer
- .1 manometer 0-16 bar for the visual control of supplied pressure
- .1 emergency supplying inlet as per NIST standard
- .1 safety valve set at 12 bar



