





## SECOND STAGE PRESSURE REDUCER

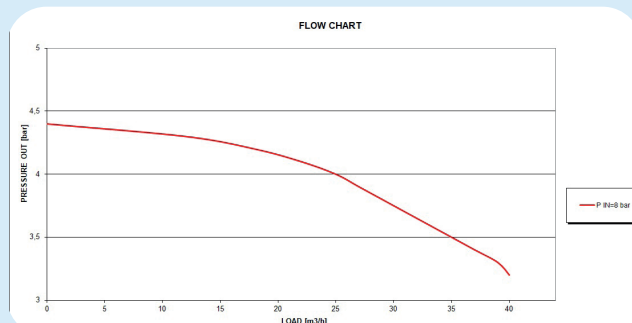

**The second stage pressure line reducer for medical gases are conceived and realized according to MDD 93/42/EEC and related ISO EN 7396-1 and ISO EN 10524-1 norms**

The second stage reducer is conceived and manufactured for regulating the pressure from primary network to set pressure of terminal units (4-5 bar as per EN ISO 7396-1 or 8-9 bar for Nitrogen and Instrument Air)

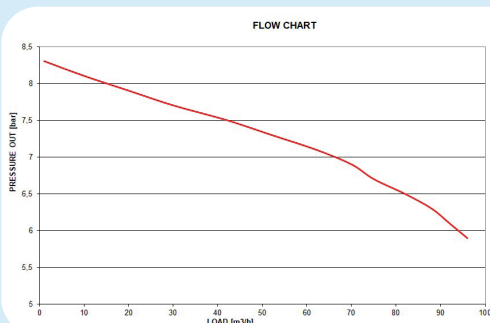

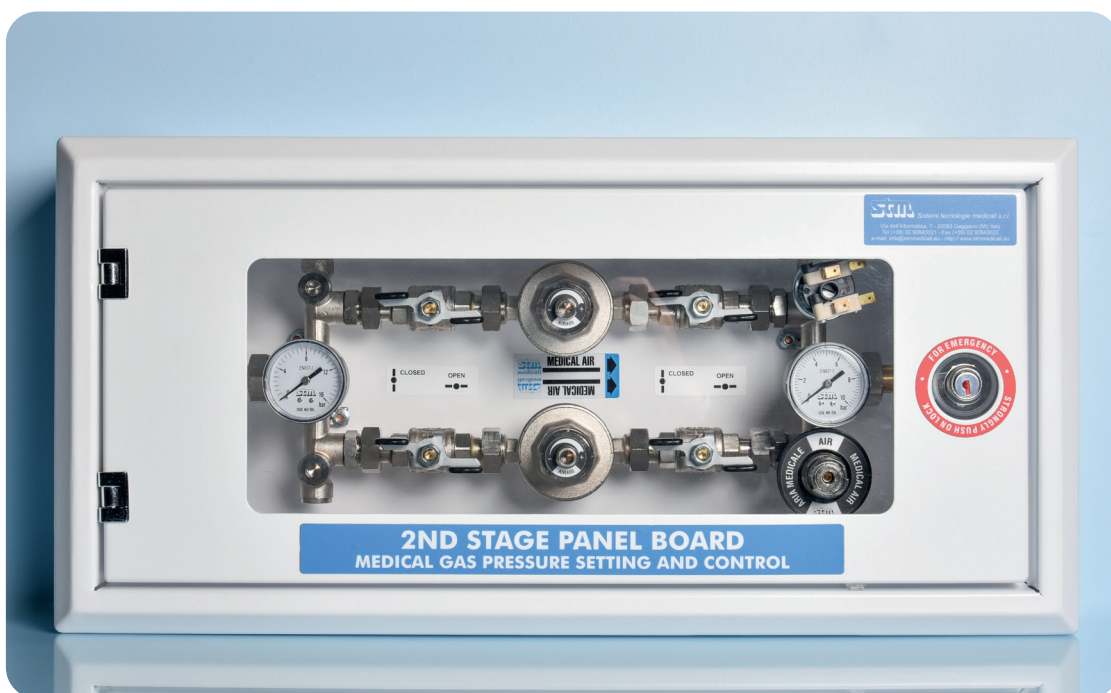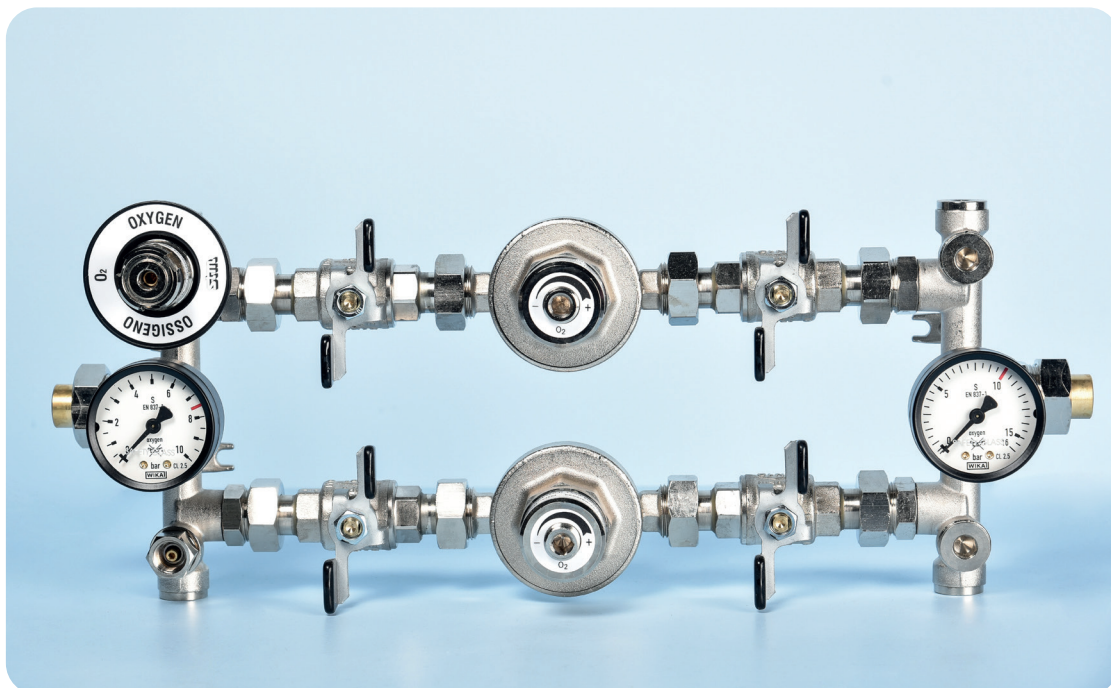

**Double second stage reducer in by-pass solution**  
**In/out ball valve for sectioning and maintenance**  
**Manometer scale 0-16 bar and 0-10 bar**  
**Pressure switches (standard) or pressure transducer (optional)**



STM SISTEMI TECNOLOGIE MEDICALI S.r.l. - VIA DELL'INFORMATICA, 7 - 20083 VIGANO DI GAGGIANO (MI) Italy  
 Tel: +39 02/90843021 - FAX: +39 02/90843022 - info@stmmedicali.eu - www.stmmmedicali.eu